


# BARI-FOAM™ CUSHIONS

Contoured bottom for sling-seat wheelchairs

Reorder #	Size(Inches)	Reorder #	Size(Inches)
☐ 754810	20 x 16 x 4	☐ 754821	24 x 20 x 4
☐ 754812	20 x 18 x 4	☐ 754823	26 x 18 x 4
☐ 754814	22 x 16 x 4	☐ 754825	26 x 20 x 4
☐ 754816	22 x 18 x 4	☐ 754827	28 x 18 x 4
☐ 754818	24 x 18 x 4	☐ 754829	28 x 20 x 4

Custom sizes are available. Please call **Skil-Care at 800-431-2972.**

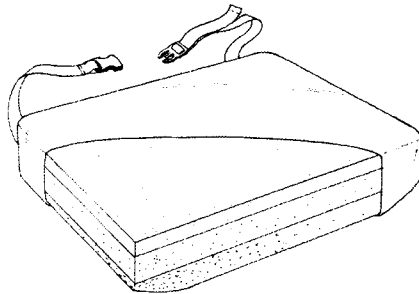

## Special construction

Skil-Care's Bari-Foam™ Cushions provide maximum comfort and pressure reduction for bariatric wheelchair users. The three-part design incorporates a top layer of temperature sensitive visco-elastic foam that molds to the user to eliminate pressure "hot spots." The center layer of high-density foam offers comfort, cushioning, and protection against bottoming-out. The firm base distributes the user's weight for efficient pressure reduction. It also eliminates sling-seat hammocking to improve sitting posture and reduce pressure on the spine and coccyx.

## Low-Shear II cover fabric

This liquid-resistant, double-stretch, breathable polyurethane fabric with laminated polyester backing reduces skin-damaging shear and friction.

## Important safety notice

Be certain to secure the cushion's safety strap to the wheelchair frame to prevent the cushion from sliding off the seat.

## Easy care

The Low-Shear II cover wipes clean with a damp cloth. Mild spray-type cleansers and disinfectants may be used. **Do not launder the cover.**