


BARI-FOAM™ CUSHIONS

Flat bottom for solid-seat wheelchairs

Reorder #	Size(Inches)	Reorder #	Size(Inches)
<input type="checkbox"/> 754833	20 x 16 x 4	<input type="checkbox"/> 754843	24 x 20 x 4
<input type="checkbox"/> 754835	20 x 18 x 4	<input type="checkbox"/> 754845	26 x 18 x 4
<input type="checkbox"/> 754837	22 x 16 x 4	<input type="checkbox"/> 754847	26 x 20 x 4
<input type="checkbox"/> 754839	22 x 18 x 4	<input type="checkbox"/> 754849	28 x 18 x 4
<input type="checkbox"/> 754841	24 x 18 x 4	<input type="checkbox"/> 754851	28 x 20 x 4

Custom sizes are available. Please call Skil-Care at 800-431-2972.

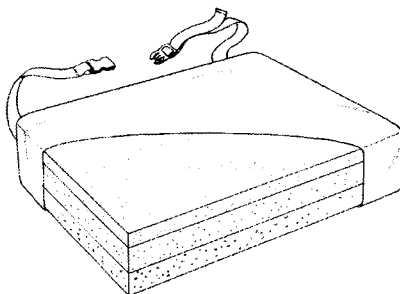

Special construction

Skil-Care's Bari-Foam™ Cushions provide maximum comfort and pressure reduction for bariatric wheelchair users. The three-part design incorporates a top layer of temperature sensitive visco-elastic foam that molds to the user to eliminate pressure "hot spots." The center layer of high-density foam offers comfort, cushioning, and protection against bottoming-out. The firm base distributes the user's weight for efficient pressure reduction.

Low-Shear II cover fabric

This liquid-resistant, double-stretch, breathable polyurethane fabric with laminated polyester backing reduces skin-damaging shear and friction.

Important safety notice

Be certain to secure the cushion's safety strap to the wheelchair frame to prevent the cushion from sliding off the seat.

Easy care

The Low-Shear II cover wipes clean with a damp cloth. Mild spray-type cleansers and disinfectants may be used. **Do not launder the cover.**