


## ECONO-GEL PAD

Reorder #	Cover	Size
751150	Vinyl w/Polyester	16" x 16" x 2"
751151	Vinyl w/Polyester	18" x 16" x 2"
751154	Vinyl Only	16" x 16" x 2"
751155	Vinyl Only	18" x 16" x 2"

- **Polyester outer cover** reduces skin shear. Two-tone fabric serves as a reminder to reposition patient. Removes for laundering.

- **Fluid-proof vinyl cover** protects against incontinence. Wipes clean for easy care.

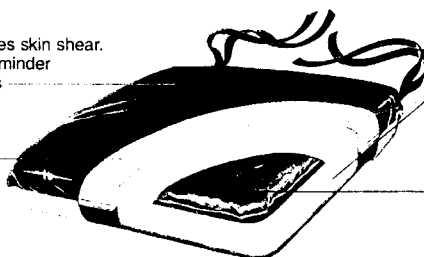

- **Safety ties** keep the cushion in place.

- **High-density foam** layers on top and bottom increase comfort and stability.

- **Gel-filled pod** between foam layers provides pressure-relieving flotation.

### Important Safety Reminder

Each Econo-Gel Pad has safety ties that must be secured to the wheelchair frame to prevent the cushion from slipping off the seat.

### Easy Care

The vinyl cover may be wiped clean using a damp cloth and spray-type cleanser. A mild disinfectant may be used. The polyester cover can be laundered and dried at temperatures that do not exceed 180°F (82°C). Do not use bleach.