


Perimeter-Guard Gel-Foam Mattress

With low-pressure gel coccyx-support pad

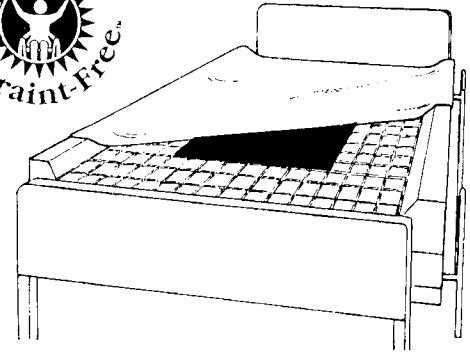

Reorder #	Size	Reorder #	Size
<input type="checkbox"/> 558055	6" x 36" x 72"	<input type="checkbox"/> 558057	6" x 36" x 80"
<input type="checkbox"/> 558056	6" x 36" x 76"	<input type="checkbox"/> 558058	6" x 36" x 84"

Specially designed for nursing homes, hospitals, and homecare. Includes these clinically tested features . . .

Perimeter-Guard Bolsters. Raised, high-density foam bolsters along the sides of the mattress offer the margin of safety patients need when bed rails cannot be used.

Pressure-Check™ Foam. Nearly 400 separate foam sections™ react individually to achieve the weight distribution patients require for reduced interface pressure.

Ultra-Soft Heel Zone. Lowers pressure for sensitive heels.


High-Density Support Foam. Provides uniform, compression-resistant support for patient comfort and extended mattress life.

Low-Shear Cover. Supple, stretch-knit fabric reduces friction to protect against skin shearing. Backed with incontinence-proof barrier. Bacterio-static and completely launderable for maximum infection control. (See laundering instruction tag on cover.)

Zipper Flap. Protects foam against moisture and contaminants.

Low-Pressure Gel Coccyx-Support Pad.

Substantially lowers interface pressure beneath the coccyx. The pad contains a waterbased gel that moves with the patient to decrease skin shearing. Its cooling effect reduces skin-damaging perspiration.

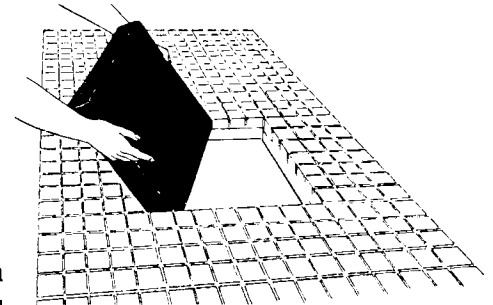

All Perimeter-Guard Mattress components conform to California Bulletin 117 flammability standards.

Bonus Values

(1) Use the gel pad as a wheelchair flotation cushion. Simply open the mattress cover, remove the gel pad, and place in on the wheelchair seat.

(2) The Perimeter-Guard Gel-Foam Mattress is packaged with a foam insert. This may be used as substitute for the gel pad for patients requiring less pressure relief.

