


SLIDE-GUARD WEDGE


Re-Order #	Description	Foundation	Cover
<input type="checkbox"/> 751310	Gel-Foam	Firm	Polyester
<input type="checkbox"/> 751320	Gel-Foam	Soft	Polyester
<input type="checkbox"/> 751350	Foam	Firm	Polyester
<input type="checkbox"/> 751370	Foam	Soft	Polyester
<input type="checkbox"/>			

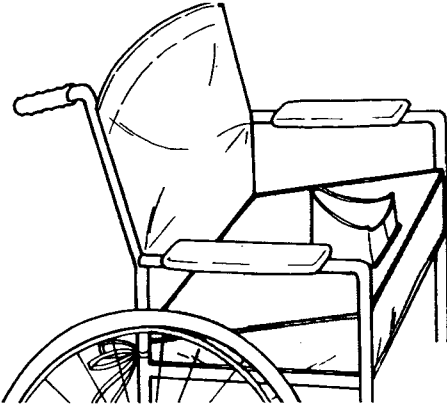

PURPOSE

Skil-Care's SLIDE-GUARD WEDGE CUSHION provides double slider control for wheelchair patients. This wedge-shaped cushion lifts the patient's legs and shifts his/her weight toward the back of the seat, reducing the tendency to slip forward. The pommel offers added protection against sliding. In addition, the pommel acts as a leg abductor.

Each Slide-Guard Wedge has ties for the purpose of securing the cushion to the wheelchair. For patient safety, these ties must be fastened to the chair.

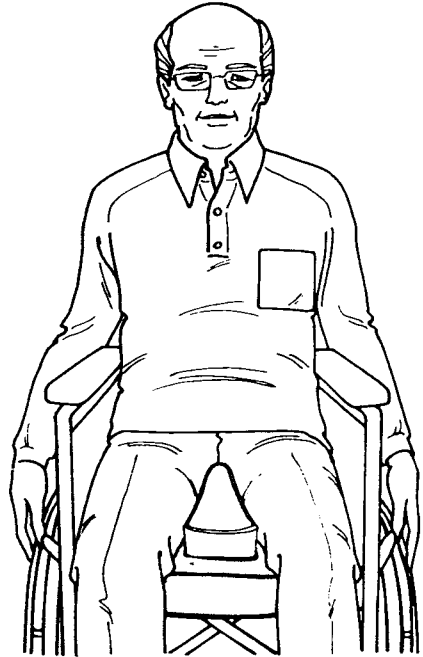

COVER LAUNDERING AND CUSHION CARE

The polyester covers remove for laundering. Do not exceed 180°F (82°C) when washing or drying these covers. Bleaching is not recommended. The vinyl cushion may be cleaned with a mild spray-type cleaner and damp cloth or sponge.