


VISCO-FOAM CUSHIONS

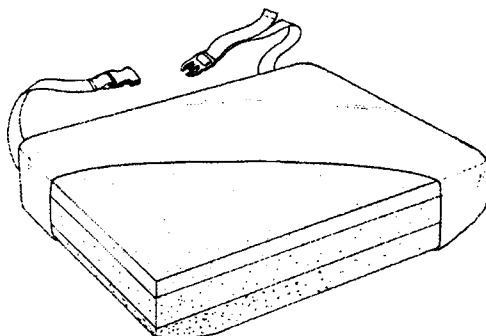

Reorder #	Size
<input type="checkbox"/> 754750	16" x 16" x 4"
<input type="checkbox"/> 754755	18" x 16" x 4"

Special construction for pressure relief and comfort

Skil-Care's Visco-Foam Cushion provides pressure relief and comfort through a unique three-part design. The cushion features a top layer of temperature-sensitive foam that molds to the user to eliminate pressure "hot spot." The center layer of high-density foam offers comfort, cushioning, and protection against bottoming out. The firm base distributes the user's weight for efficient reduction. It also eliminates sling-seat hammocking to improve sitting posture and reduce pressure on the spine and coccyx.

Low-Shear II cover fabric

Be certain to secure the cushion's safety strap to the wheelchair frame in order to prevent the cushion from sliding off the seat.

Easy Care

The Low-Shear II cover wipes clean with a damp cloth. Mild spray-type cleansers and disinfectants may be used. **Do not launder this cover.**